Tubes Radiogènes à Anode Tournante  
Drehanoden - Röntgenröhre  
Tubos de Rayos-X con Ánodo Giratorio

Note: Document originally drafted in the English language.

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**Note:**  
1. Limits are based on maximum track rating except for the following codes:  
a - Limited by available heat storage.  
b - Limited by window heating.  
c - Limited by filament emission.  
2. H.S. = Heat Storage  
kV = Tube Voltage

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**Note:**  
Rating charts reflect maximum tube performance. Tube operation is ultimately limited by system software.

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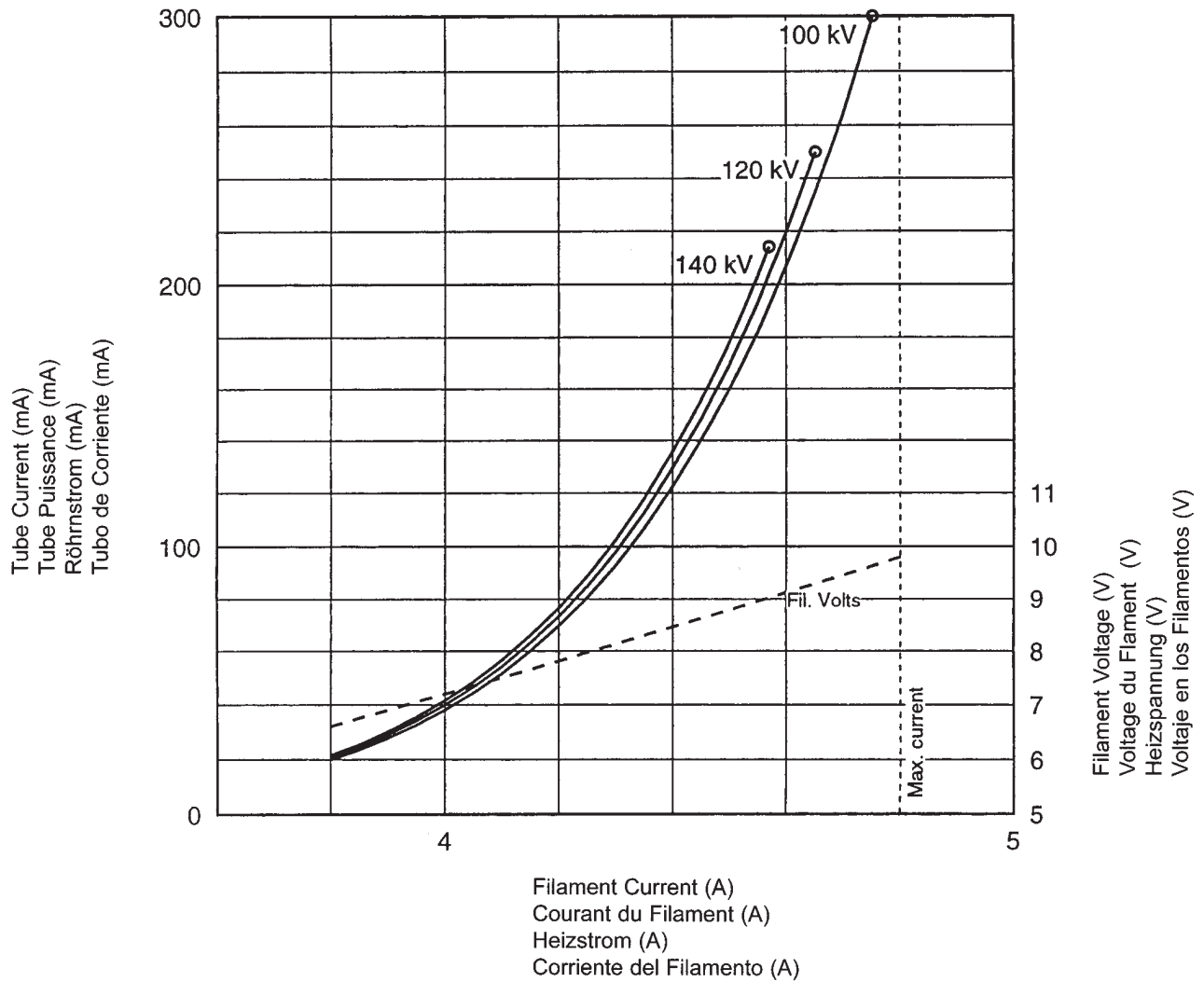
Le Gaine B-240H

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Encaje de B-240H

Maximum Peak Voltage .....	150 kV
Anode to Ground .....	75 kV
Cathode to Ground .....	75 kV
Maximum X-ray Tube Assembly Heat Content .....	1.5 MJ (2.0 MHU)
Maximum Continuous Heat Dissipation (Includes stator heat) .....	3.7 kW (5.18 kHU/sec)
Maximum Heat Exchanger Dissipation .....	5.0 kW (7.0 kHU/sec)
Focal Point Position (Central Ray) Within 1mm (X, Y Direction from the center of radiation port.)	
X-Ray Tube Assembly Permanent filtration .....	1.0 mm Al IEC 60522
Loading Factors for Leakage Radiation .....	150 kV, 23 mA
Federal Standard High Voltage Cable .....	72
Ambient Air Temperature Limits for Operation .....	5°C to 40°C
Temperature Limits for Storage and Transport .....	-20°C to +70°C
Humidity .....	10% - 90%
Weight: Housing .....	33.1 kg (73 lbs)
Heat Exchanger .....	15.9 kg (35 lbs)
IEC Classification .....	Class 1
Safety Devices: Thermal Switch Normally Closed Contact .....	Opening at 85°C
Flow Switch - Normally Open contact Contacts close with adequate oil flow.	
Filament Frequency Limits .....	50 HZ - 25 kHz
Power Supply .....	DC

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
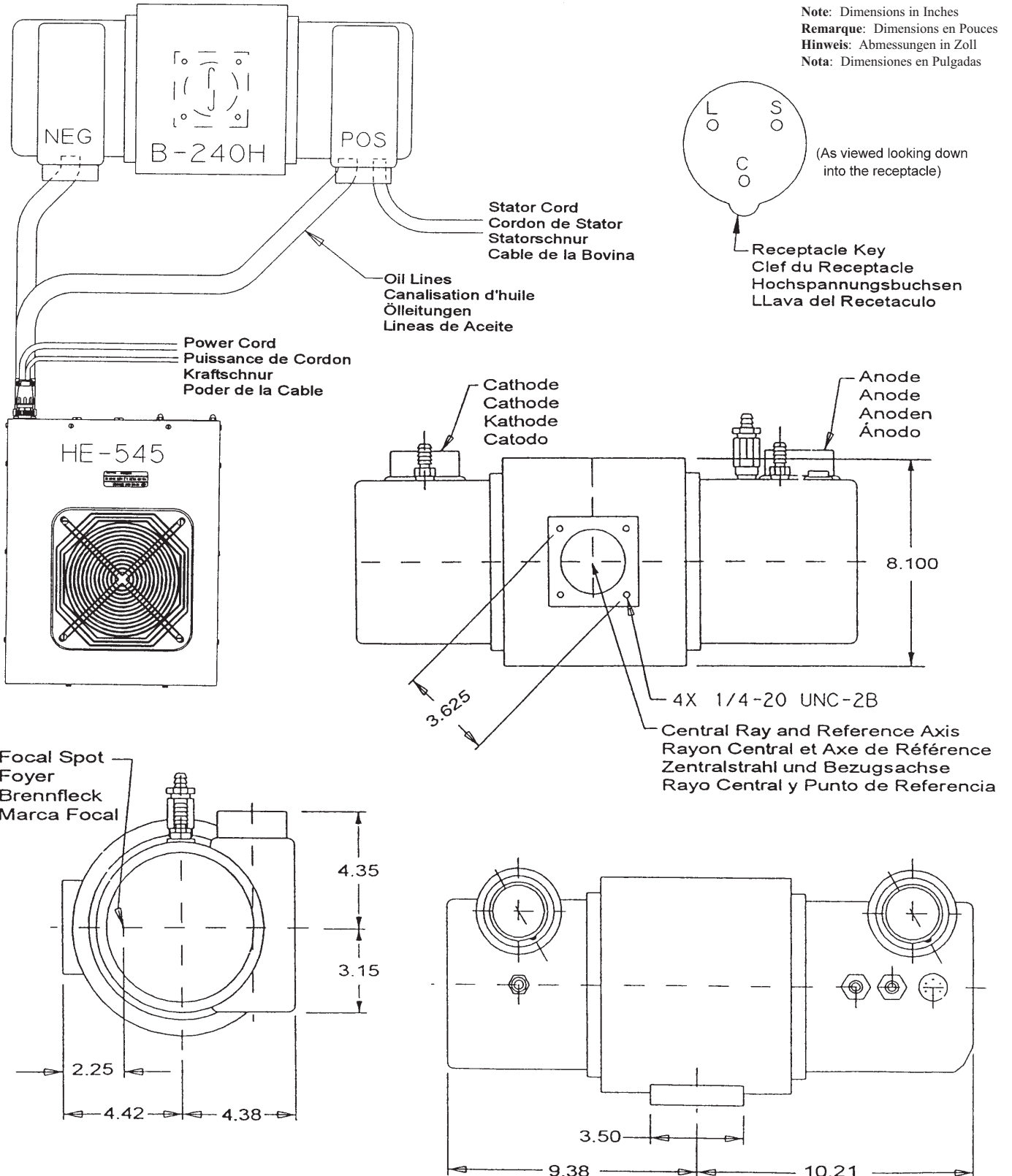
**TUBE HOUSING ASSEMBLY HEATING AND COOLING**

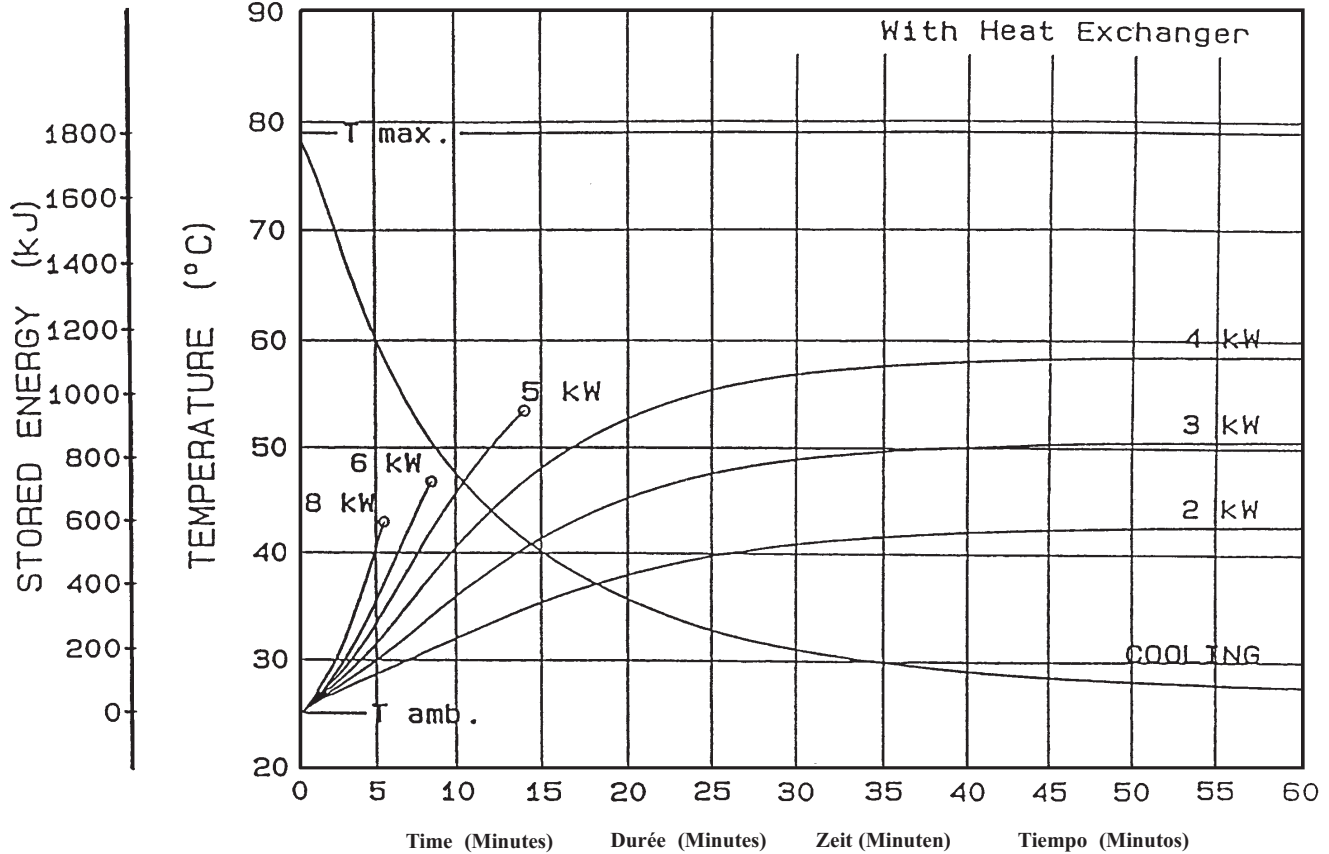**Note:**  
1. Heat inputs into housing include tube power, filament power, and stator power.  
2. Heating curves based on no restrictions of natural convection around tube housing assembly.  
3. Heating and cooling curves reflect maximum tube performance. Tube operation is ultimately limited by system software control.

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
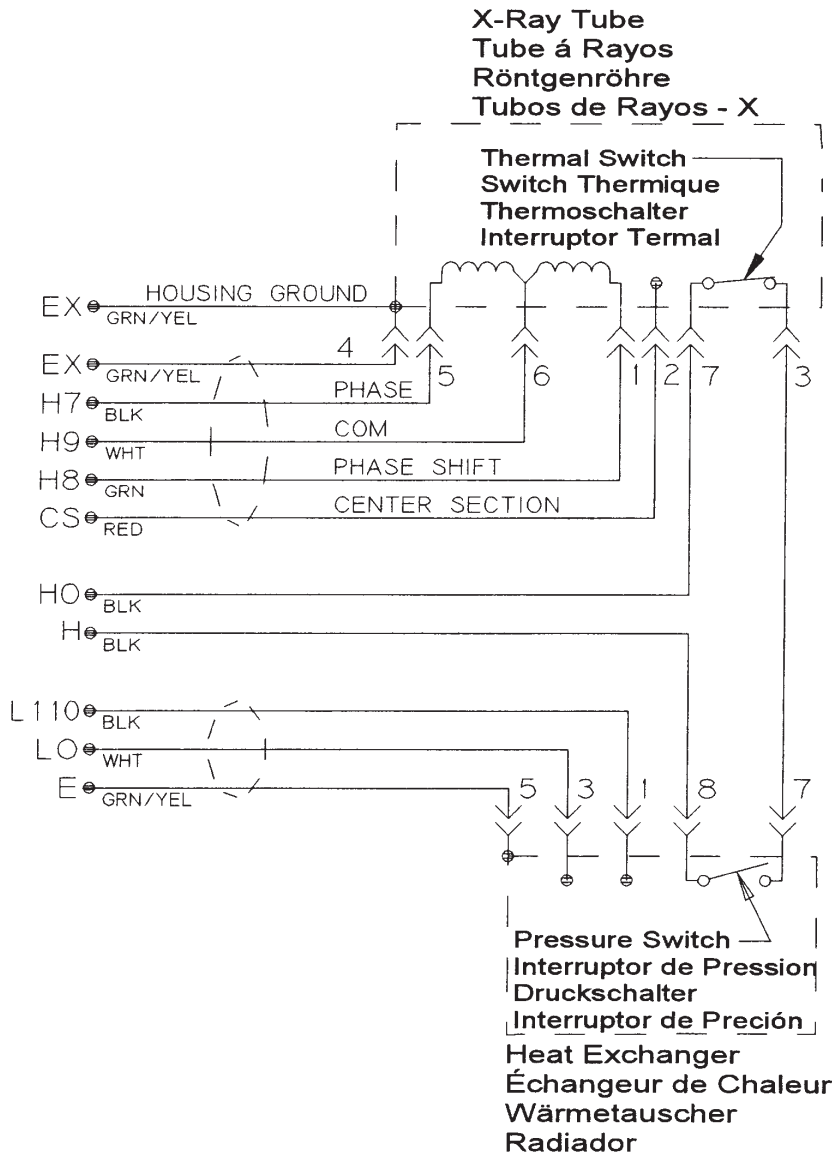
<b>Stator Type: "R"</b>	
<b>Stator Coil Resistance:</b>	
Black to White	14 Ohms ±15%
Green to White	46 Ohms ±15%
<b>Starter Voltage:</b>	
50/60 Hz	<b>Start</b> 200 VAC <b>Run</b> 80 VAC
<b>Time to Full Speed:</b>	
50/60 Hz	0 - 2700 RPM 9.6 Sec.
<b>X-Ray Tube Assembly:</b>	
GS-3072/B-240H	IEC 60601-2-28

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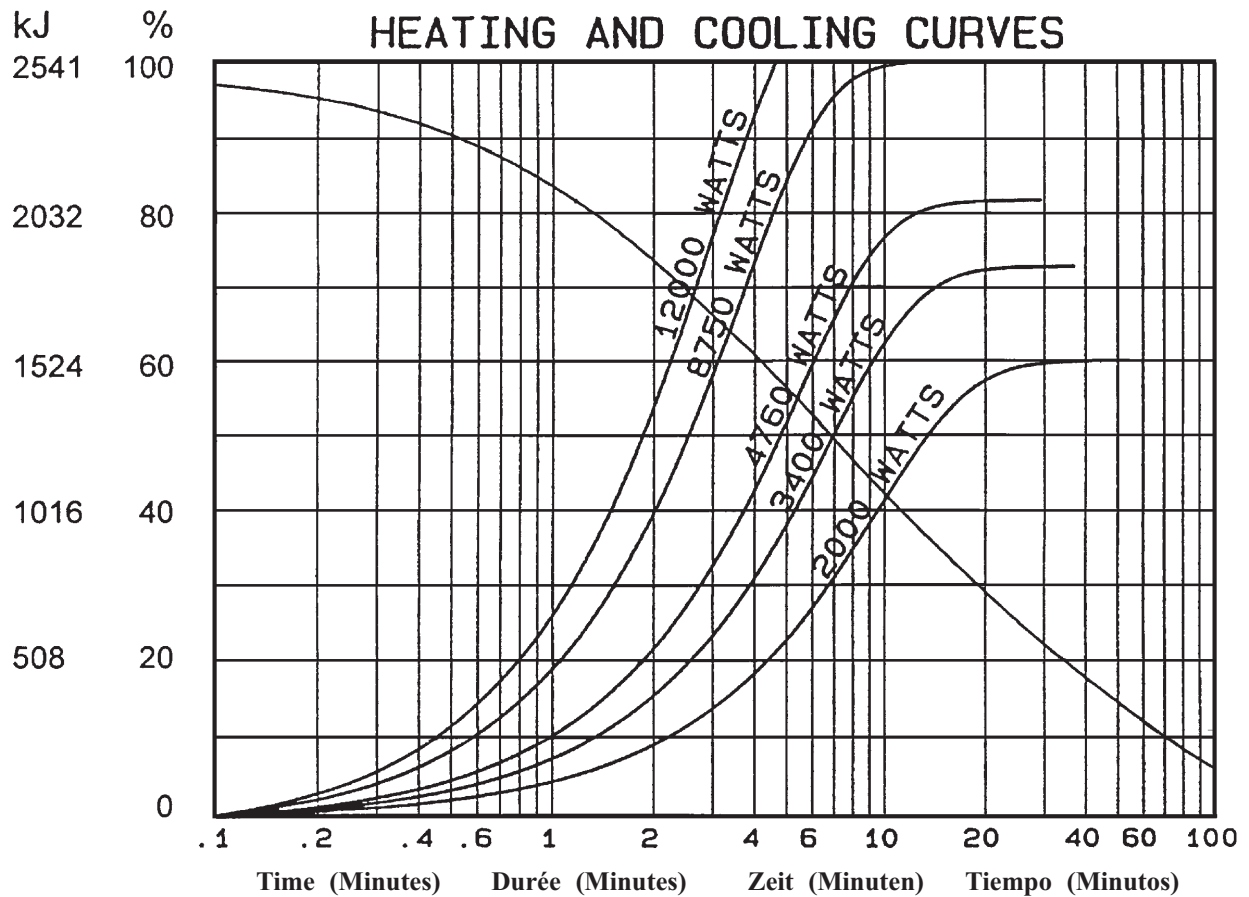**Note:**  
1. Heating and cooling curves reflect maximum tube performance. Tube operation is ultimately limited by system software control.

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